About Signia

Signia is one of the leading hearing aid brands in the world. Our goal is to enhance human performance through iconic innovations that shape the prescription audio category. We refuse to accept hearing loss as a limitation. It's not about correcting a loss, it's about gaining an edge. Not just performing, performing brilliantly. Not just hearing again, but hearing better than ever. Hearing every word, beat, and breath. With Signia, you can Be Brilliant.

Find out more on signia-pro.com

Or ask your hearing care professional for more information on Signia hearing aids, useful accessories and the Signia app.

Download on the App Store | GET IT ON Google Play

Compatible with Android and iOS. The app is free of charge.

The information in this document contains general descriptions of the technical options available, which do not always have to be present in individual cases and are subject to change without prior notice.

Apple, the Apple Logo, App Store, and iPhone are trademarks of Apple Inc. Android is a trademark of Google LLC. The Bluetooth® word mark and logos are owned by the Bluetooth SIG, Inc. and any use of such marks by WSAUD A/S is under license.The 'Qi' symbol is a trademark of the Wireless Power Consortium. Any other trademarks and trade names are those of their respective owners.

Find us on

- signia-pro.com
- facebook.com/signiahearing
- youtube.com/signiahearing
- instagram.com/signiahearing

Order No. 05681-99C1-7600
© 08.2025 WSAUD A/S



signia



Be Brilliant™

Signia is one of the leading hearing aid brands in the world. Our goal is to enhance human performance through iconic innovations that shape the prescription audio category. We refuse to accept hearing loss as a limitation. It's not about correcting a loss, it's about gaining an edge. Not just performing, performing brilliantly. Not just hearing again, but hearing better than ever. Hearing every word, beat, and breath.

With Signia, you can Be Brilliant.

Contents

Signia Integrated Xperience	4
Styletto IX	8
Pure Charge&Go IX family	9
Active Pro IX /Active IX	10
Motion Charge&Go IX	11
Silk Charge&Go IX	14
Insio IX family	15
IX Feature Overviews	18
Signia Augmented Xperience	24
Styletto AX	26
Pure Charge&Go AX / Pure Charge&Go T AX	27
Insio Charge&Go AX	28
Pure 312 AX	29
AX Feature Overview	30
Signia Xperience	32
Xperience Feature Overview	32
Motion Charge&Go X, P X, SP X	33
Silk X	36
Pure 312 X	37
Sirion Connect	38
Orion Charge&Go	42
Intuis 4	46
Intuis 3	48
Basic	52
Color overview	58
EarWear	60
Chargers	64
Accessories	66
Signia app	68
Fitting	70

Signia IX

Unleash the power of conversation

Performance levels

- 7IX
- 5IX
- 3IX
- 2IX
- 1IX

Meet the world's first hearing aid with a multi-stream architecture capable of accurately pinpointing multiple moving speakers in real-time conversation.

IX Integrated Xperience

Improved hearing is about more than just sound amplification. It is about empowering individuals to immerse themselves in their surroundings, create meaningful connections, and cherish every moment of life's experiences.

Leveraging our profound understanding of dynamic conversations, including the head movements of the wearer and changing positions of speakers, we crafted Signia Integrated Xperience. This groundbreaking technology ensures uninterrupted conversation engagement in noise, addressing a critical challenge many hearing aid wearers face.

Our commitment is to support wearers in leading active lives, unshackled from the limitations of traditional hearing aids and free to contribute and foster deeper connections.

Experience the transformative power of real-time, dynamic conversation enhancement with Signia Integrated Xperience.

SLIM-RIC



Styletto IX

RIC



Pure Charge&Go IX / TIX / BCT IX

BTE



Motion Charge&Go M IX

ITE



Active Pro IX



Insio Charge&Go CIC IX



Silk Charge&Go IX



Insio IIC/CIC IX



Insio ITC IX



Insio ITE IX



RealTime Conversation Enhancement

The science of keeping up with conversations

Signia Integrated Xperience's revolutionary RealTime Conversation Enhancement analyzes, augments, and adapts to the dynamic flow of conversations. This lets wearers enjoy a seamless and enriched hearing experience.

Analyze
Processes 192,000 data points every second to identify speakers and analyze conversation dynamics, authentically reflecting real-life conversations.

Augment
Uses data gathered to control multiple focus streams that lock on to and dynamically highlight each speaker's voice amidst the noise. Meanwhile, Augmented Focus technology[2] processes background sound to ensure a comfortable, immersive conversational experience without distractions.

Adapt
Our patented multi-stream architecture updates 1000 times per second, merging multiple focus streams into a dynamic auditory space that adapts to evolving conversations. This allows for real-time tracking and enhancement of multiple moving speakers, ensuring wearers can fully engage in group discussions.



A conversation starter



Key features

- RealTime Conversation Enhancement
- Augmented Focus
- Own Voice Processing 2.0
- Signia Assistant
- MFi, LE Audio, ASHA
- Convenient rechargeability
- Signia app
- CROS compatibility
- Sleek, stylish design

Performance levels

7IX
5IX
3IX
2IX
1IX

sDemo

Styletto IX

Sleek, sophisticated, and available in a range of distinctive colors, Styletto IX hearing aids blend superior IX technology with style. Enjoy great hearing and a great look – and the added convenience of recharging on the go.

Also available as CROS Styletto IX.



Charging options
The sleek, portable Styletto Charger has an integrated power bank that provides up to five additional charges. The Styletto Dry&Clean charger not only charges the hearing aids but also offers drying and UV sanitizing technology. See pages 60-61 for more information on chargers.

SLIM-RIC

Model			
Receiver type	miniReceiver 3.0 S	miniReceiver 3.0 M	miniReceiver 3.0 P
			
Performance			
Matrix (2 ccm)	46/110 dB SPL	60/119 dB SPL	65/122 dB SPL
Battery			
Type	Lithium-ion power cell	Lithium-ion power cell	Lithium-ion power cell

Catch all sides of the conversation



Pure Charge&Go IX



Pure Charge&Go T IX



Pure Charge&Go BCT IX

Key features

- RealTime Conversation Enhancement
- Augmented Focus
- Own Voice Processing 2.0
- Signia Assistant
- Option for Bluetooth Classic or MFi, LE Audio, ASHA
- Convenient rechargeability
- Signia app
- CROS compatibility
- Telecoil options available

Performance levels

7IX
5IX
3IX
2IX
1IX

sDemo

Pure Charge&Go IX family

With the addition of Pure Charge&Go BCT IX to the Pure Charge&Go IX family, you can now offer all your clients the full power of IX regardless of Bluetooth preference.

Pure Charge&Go BCT IX offers Bluetooth Classic connectivity, an integrated telecoil, and an LED.

Pure Charge&Go IX / T IX offer connectivity via MFi, ASHA, and LE Audio, and Pure Charge&Go T IX also has an integrated telecoil.

All three options offer plenty of runtime as well as CROS compatibility.



Charging options
The Pure Portable Charger provides convenient charging on the go. The Pure Dry&Clean Charger also offers UV cleaning and drying to keep the hearing aids in peak condition. And the small, reliable standard Pure Charger meets basic desktop charging needs. Also compatible with the Multi Charger. See pages 60-61 for more information on chargers.

RIC

Model				
Receiver type	miniReceiver 3.0 S	miniReceiver 3.0 M	miniReceiver 3.0 P	miniReceiver 3.0 HP
				
Performance				
Matrix (2 ccm)	46/110 dB SPL	60/119 dB SPL	65/122 dB SPL	75/131 dB SPL
Battery				
Type	Lithium-ion power cell	Lithium-ion power cell	Lithium-ion power cell	Lithium-ion power cell

Reconnect with the conversation



Active Pro IX / Active IX

Welcome to a new era of hearing with Signia Active Pro IX – an earbud on the outside, with Signia's game-changing prescription-grade technology on the inside - including RealTime Conversation Enhancement. Fits instantly, enables hands-free calls and streaming, and resists sweat and dust.

Also available as a CROS Active IX.

Performance levels

- Active IX Pro
- Active IX
- sDemo

Key features

- RealTime Conversation Enhancement
- Augmented Focus
- Signia Assistant
- Android & iPhone connectivity
- Convenient rechargeability
- Signia app
- CROS compatibility
- Instant fit earbud



Active IX Charger
Signia Active Pro IX / Active IX comes in a pocket-size charging case with an integrated battery bank that can charge a pair of devices three times before being recharged. The charging case can be charged wirelessly or via USB-C cable.

	Instant Fit ITE
Model	
	
Performance	
Matrix (2 ccm)	47/117 dB SPL
Battery	
Type	Lithium-ion power cell

Bring everyone into the conversation



Motion Charge&Go M IX

With the Motion Charge&Go IX BTE family, we introduce the broadest BTE range with proven technology to enhance speech understanding in noise - RealTime Conversation Enhancement.

Performance levels

- 7IX
- 5IX
- 3IX
- sDemo

Key features

- RealTime Conversation Enhancement
- Augmented Focus
- OVP Own Voice Processing 2.0
- Signia Assistant
- MFi, LE Audio, ASHA
- Convenient rechargeability
- Signia app
- CROS compatibility
- Telecoil



Convenient rechargeability
Compatible with the Desktop Charger Motion M, a small desktop charger with LEDs to indicate charging status, as well as with the Multi Charger.

	BTE		
Model			
	Earhook	ThinTube 3.0	ThinTube 3.0 P
			
Performance			
Matrix (2 ccm)	60/133 dB SPL	60/125 dB SPL	63/126 dB SPL
Battery			
Type	Lithium-ion power cell	Lithium-ion power cell	Lithium-ion power cell

Connect with conversation



Motion Charge&Go P IX

Performance levels

- 7IX
- 5IX
- 3IX
- sDemo

The Motion Charge&Go BTE IX family embodies our belief that everyone should be able to contribute to the conversation. That's why RealTime Conversation Enhancement, our unique technology for enhancing conversations in noise, is available in BTEs for all levels of hearing enhancement.

Key features

- RealTime Conversation Enhancement
- Augmented Focus
- Own Voice Processing 2.0
- Signia Assistant
- MFi, LE Audio, ASHA
- Convenient rechargeability
- Signia app
- CROS compatibility
- Telecoil



Convenient rechargeability

Charging options include the Desktop Charger Motion P, the Dry&Clean Motion P which also dries and sanitizes the hearing aids, and the Multi Charger.

	BTE		
Model			
	Earhook	ThinTube 3.0	ThinTube 3.0 P
			
Performance			
Matrix (2 ccm)	77/135 dB SPL	66/130 dB SPL	70/131 dB SPL
Battery			
Type	Lithium-ion power cell	Lithium-ion power cell	Lithium-ion power cell

Superpower your hearing



Motion Charge&Go SP IX

Performance levels

- 7IX
- 5IX
- 3IX
- sDemo

All your clients - even those requiring the most amplification - can benefit from our groundbreaking technology for enhancing conversations, RealTime Conversation Enhancement. And they'll never miss a moment with the longest battery runtime of any superpower device, up to 65 hours from a single charge.

Key features

- RealTime Conversation Enhancement
- Augmented Focus
- Own Voice Processing 2.0
- Signia Assistant
- MFi, LE Audio, ASHA
- Convenient rechargeability
- Signia app
- CROS compatibility
- Telecoil



Convenient rechargeability

Charging options include the Desktop Charger Motion SP, the Dry&Clean Motion SP which also dries and sanitizes the hearing aids, and the Multi Charger.

	BTE
Model	
	Earhook
	
Performance	
Matrix (2 ccm)	82/140 dB SPL
Battery	
Type	Lithium-ion power cell

The only thing that changes is your hearing



Performance levels

7IX
5IX
3IX
2IX
1IX

sDemo


Silk Charge&Go IX

Comfortable and discreet, this tiny CIC hearing aid provides the convenience of on-the-go rechargeability. It is packed with advanced sound features to support you in any environment including its unique Binaural OneMic Directionality, which offers unmatched conversation support in noisy environments for a hearing aid of its size.

Also available as a CROS Silk Charge&Go IX.

Key features

- Binaural OneMic Directionality
- On-the-go rechargeability
- Instant fit
- Discreet CIC fit
- Signia app
- CROS compatibility

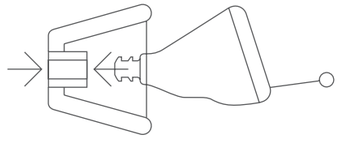

Instant Fit
EarWear 3.0 sleeves come in four different sizes, with closed and open versions to suit most ear anatomies and hearing needs. Also compatible with custom Earmold 3.0.



Silk IX Charger
Pocket-sized charging case with integrated power bank that can provide up to 3 full charges before requiring a recharge. Can be charged wirelessly.

Model	Instant Fit ITE
	
Performance	
Matrix (2 ccm)	50/114 dB SPL
Battery	
Type	Lithium-ion power cell

Invisibly Powerful



Key features: one-microphone models

- Binaural OneMic Directionality
- On-the-go rechargeability
- Custom fit
- Discreet CIC fit
- Signia app
- CROS compatibility

Performance levels

7IX
5IX
3IX

sDemo

Insio Charge&Go CIC IX

Thanks to wide bandwidth, high speed e2e 4.0, Insio Charge&Go CIC IX is uniquely able to create a directional beam in noise, ensuring the wearer hears the person in front of them without distractions. Our most discreet rechargeable custom devices yet, they provide up to 35 hours of battery runtime, and a pocket-sized charging case provides the option to charge on the go.



Portable Charger Insio CIC
Pocket-sized charging case with integrated power bank that can provide up to three full charges for a pair of hearing aids. Can be charged wirelessly.

Model	CIC		
			
Performance			
Matrix (2 ccm)	50/113 dB SPL	55/118 dB SPL	60/124 dB SPL
Battery			
Type	Lithium-ion power cell	Lithium-ion power cell	Lithium-ion power cell

Engineered for conversations

Key features: two-microphone models

- RealTime Conversation Enhancement
- Augmented Focus
- Custom fit
- Signia app
- CROS compatibility



Performance levels

7IX
5IX
3IX
2IX
1IX

Insio ITC/ITE IX

Offering the option for two microphones, Insio IX ITC* and ITE deliver the full benefit of Signia's groundbreaking RealTime Conversation Enhancement technology in a comfortable custom design.

Custom ITE

Model	ITC S	ITC M /ITE M	ITC P / ITE P
			
Performance			
Matrix (2 ccm)	50/113 dB SPL	55/118 dB SPL	60/124 dB SPL
Battery			
Size	13	13	13
Power consumption	1.0 mA**	1.1 mA**	1.1 mA**

*) Insio IX ITC is also available in a one-microphone version, which provides Binaural OneMic Directionality instead of RealTime Conversation Enhancement.
**) Max. values

Engineered for conversations

Key features: one-microphone models

- Binaural OneMic Directionality
- Custom fit
- Discreet IIC or CIC fit
- Signia app
- CROS compatibility



Performance levels

7IX
5IX
3IX
2IX
1IX

Insio IIC/CIC IX

With a new design, and built on our IX architecture, our popular Insio IICs and CICs are better than ever. Insio IX is the only custom one-mic hearing aid offering binaural directionality for focused conversation, even in noisy environments.

Custom ITE

Model	IIC / CIC S	CIC M	CIC P	
	I3* 	I3* 	I3* 	T5** 
Performance				
Matrix (2 ccm)	50/113 dB SPL	55/118 dB SPL	60/124 dB SPL	
Battery				
Size	10	10	10	
Power consumption	1.0 mA***	1.0 mA***	1.0 mA***	

*) I3 faceplate available for IIC and CIC in e2e wireless and non-wireless versions. Non-wireless versions programmed with I3 Adapter, wireless can be programmed with ConnexxAir.
**) T5 faceplate available if push-button is desired. Non-wireless versions programmed with Flex Connector. Wireless versions can be programmed with ConnexxAir.
***) Max. values

Feature overview

Pure Charge&Go IX family, Styletto IX, Active Pro IX, Motion Charge&Go IX family†, Insio ITC††/ITE IX

Features	Description	Performance levels: Styletto IX, Pure Charge&Go IX family, Motion Charge&Go IX family†, Insio IX ITC††/ITE IX					Performance levels: Active Pro IX	
		7 IX	5 IX	3 IX	2 IX	1 IX	Active Pro	Active
Conversation management								
RealTime Conversation Enhancement	Provides uninterrupted engagement in noisy group conversations by seamlessly tracking the locations of all conversation partners and enhancing their speech.[1]	●●●●●	●●●●○	●●○○○			●●●●●	
Conversation Dynamics Analyzer	Analyzes conversation dynamics including turn-taking and timing, locations of pure noise without speech, and levels of background noise.	✓					✓	
Conversation Source Analyzer	Constantly scouts the acoustic environment to identify number of talkers and their location.	✓	✓	✓			✓	
Dynamic Focus Streams	Isolate and track speech independently.	3	1				3	
Conversation Booster	Dynamically enhances individual speakers.	●●●●●	●●●●○	●●○○○			●●●●●	
Augmented Focus	Separates the focus sound (speech) from surrounding sounds, creating a strong contrast.[1]	✓	✓	✓	✓	✓	✓	✓
Focus Stream	Enhances speech signals from the front hemisphere.	1	1	1	1	1	1	1
Surrounding Stream	Processes environmental sounds to preserve environmental awareness without masking speech.	1	1	1	1	1	1	1
Dynamic Soundscape Processing 3.0	Integrates inputs for precise hearing aid optimization, including (when available): acoustic sensors, motion sensors, OVP 2.0, noise detection, SNR analysis, and RealTime Conversation Enhancement.	●●●●●	●●●●○	●●●○○	●●○○○	●●○○○	●●●●●	●●○○○
Spatial SpeechFocus	Steers directional microphone beam to the dominant speech source (left, right, or behind the wearer automatically in a car situation or in a dedicated "Strolling" program).[1]	✓					✓	
SpeechFocus	Depending on the direction of the dominant speech source, Speech Focus steers the directional beam behind the wearer, or engages a surround-sound mode for speech from either side.	✓	✓				✓	
Own Voice Processing 2.0 (OVP™ 2.0)[2]	Uses real-time recognition of the wearer's voice to deliver a natural own-voice impression.[1]	✓	✓	✓				
Motion sensor[3]	Detects whether the wearer is walking, running, or sitting still to guide how sound is optimized.	✓	✓	✓			✓	
Binaural Directionality	Provides scientifically proven enhanced speech understanding in noise.[1]	✓	✓	✓			✓	
Adaptive Directionality	Improves the signal to noise ratio in noisy environments by reducing dominant noise sources from the back and sides. The polar plot is adjusted continuously in multiple frequency channels.	✓	✓	✓	✓		✓	
Automatic Directionality	Improves the signal to noise ratio in noisy environments by automatically and seamlessly applying directionality in individual processing channels.	✓	✓	✓	✓	✓	✓	✓
TruEar 360	Reproduces benefits of the natural pinna effect by highlighting the difference between front and back sources of sound.	✓	✓	✓	✓	✓	✓	✓

1) Available for bilateral fittings with e2e wireless link
2) Not available in Insio ITC/ITE IX and Active Pro IX
3) Not available in Insio ITC/ITE IX
†) Motion Charge&Go IX BTE family only available in performance level 7IX, 5IX, 3IX
††) Two-microphone models only

Features	Description	Performance levels: Styletto IX, Pure Charge&Go IX family, Motion Charge&Go IX family†, Insio IX ITC††/ITE IX					Performance levels: Active Pro IX	
		7 IX	5 IX	3 IX	2 IX	1 IX	Active Pro	Active
SpeechStabilizer	Preserves the natural envelope of speech in any situation, including very noisy situations, which preserves acoustic cues that are important for speech understanding.	✓	✓	✓	✓	✓	✓	✓
Speech and Noise Management	Reduces listening effort using sophisticated sound processing.	✓	✓	✓	✓	✓	✓	✓
e2e wireless 4.0	Controls coupling, synchronization, CROS/BiCROS, and a range of sound features.[1]	✓	✓	✓	✓	✓	✓	✓
Listening Comfort								
Extended bandwidth[4]	12 kHz available bandwidth for enhanced processing of high-frequency speech and environmental sounds.	✓					✓	
Auto EchoShield	Automatically reduces the effect of reverberations and sudden loud noises in reverberant environments.	✓					✓	
Reverberant Room Program	Dedicated program for reverberant environments.	✓	✓				✓	
SoundSmoothing™	Reduces the annoyance of sudden loud impact sound such as when handling glassware.	●●●●●	●●●●●	●●●●●	●●●○○	●○○○○	●●●●●	●○○○○
eWindScreen™	Reduces the interference of wind noise when outdoors.	●●●●●	●●●●●	●●●○○	●○○○○	●○○○○	●●●●●	●○○○○
Extended dynamic range	117 dB 24 bit processing for enhanced digital encoding of soft and loud sounds.	✓	✓	✓	✓	✓	✓	✓
Feedback Cancellation	High speed monitoring and control of feedback in individual processing channels.	✓	✓	✓	✓	✓	✓	✓
Wearer Empowerment								
Spatial Configurator	Allows the wearer to adjust the direction of the microphone beam via the Signia App.[1]	✓	✓				✓	
My WellBeing[2]	Accessed via the Signia App, wearers can get insights about how much they are using their hearing aids, their level of social interaction, their daily step count, and their overall activity level.	✓	✓	✓			✓	
Signia Assistant[2]	Offers wearers instant 24/7 troubleshooting and customized sound adjustments via an intuitive AI-enabled chat. Provides detailed information in Connexx to address each client's fine-tuning needs more precisely.	✓	✓	✓	✓	✓	✓	✓
Signia App (iOS and Android)[3]	Easy-to-use smartphone app enables remote control adjustments, audio streaming, access to My WellBeing, and support via Signia Assistant and TeleCare.	✓	✓	✓	✓	✓	✓	✓
Adaptive Streaming Volume[2]	Speech intelligibility weighted adaptation of the streaming volume depending on the level of ambient noise.	✓	✓	✓	✓	✓	✓	✓
TeleCare[3]	Provides remote services including the TeleCare Dashboard, Remote Fine Tuning, and CareChat (convenient text messages, phone calls, and video calls).	✓	✓	✓	✓	✓	✓	✓

1) Available for bilateral fittings with e2e wireless link
2) Not available in Insio ITC/ITE IX
3) Reduced features available for Insio ITC/ITE IX
4) Not available in Motion Charge&Go SP IX
†) Motion Charge&Go IX BTE family only available in performance level 7IX, 5IX, 3IX
††) Two-microphone models only

Feature overview

Pure Charge&Go IX family, Styletto IX, Active Pro IX, Motion Charge&Go IX family†, Insio ITC††/ITE IX

Features	Description	Performance levels: Styletto IX, Pure Charge&Go IX family, Motion Charge&Go IX family†, Insio IX ITC††/ITE IX					Performance levels: Active Pro IX	
		7 IX	5 IX	3 IX	2 IX	1 IX	Active Pro	Active
Connectivity								
Direct audio streaming[2]	Bluetooth Classic connectivity available in Pure Charge&Go BCT IX. MFi, LE Audio, and ASHA connectivity available in Pure Charge&Go IX / T IX, Styletto IX, Motion Charge&Go IX BTE family and Active Pro IX.	✓	✓	✓	✓	✓	✓	✓
HandsFree[2]	Enables hands-free calling. The wearer's voice is picked up by the hearing aid microphone and delivered in high quality to the conversation partner, while the call audio is streamed directly to the wearer's hearing aids. Available for connections via Bluetooth Classic, MFi, or LE Audio.	✓	✓	✓	✓	✓	✓	✓
CallControl[2]	Calls can be answered and ended by a press of the hearing aid button or by tapping the hearing aid, depending on hearing aid model. Available for connections via Bluetooth Classic, MFi, or LE Audio.	✓	✓	✓	✓	✓	✓	✓
Telecoil	Available in Pure Charge&Go BCT IX, Pure Charge&Go T IX, the Motion Charge&Go IX BTE family, and as an option in Insio ITC/ITE IX. Enables direct connection to audio loop systems and compatible telephone signals.	✓	✓	✓	✓	✓		
StreamLine Mic compatibility[1]	Stream high quality stereo sound directly from any Bluetooth device (PCs, older Apple devices, Android devices without ASHA or LE Audio, etc.). Can also be used as a remote mic.	✓	✓	✓	✓	✓	✓	✓
TV Sound" compatibility[2]	Stream high quality stereo sound directly from TV to hearing aids via Bluetooth.	✓	✓	✓	✓	✓	✓	✓
miniPocket compatibility	Optional keychain remote control.	✓	✓	✓	✓	✓	✓	✓
Tinnitus								
Ocean Waves Therapy Signal (presets)	Therapy signals that mimic the sound of the sea to provide a soothing listening experience and reduce the impact of tinnitus.	4	4	4	4		4	
Static Therapy Signal (channels / presets)	Soft, gentle therapy signals that reduce the impact of tinnitus.	20/5	18/5	16/5	12/5		20/5	
Notch Therapy[3]	Signia's unique, patented tinnitus therapy which is proven especially effective against tonal tinnitus.	✓	✓	✓	✓		✓	

1) Not compatible with Insio ITC/ITE IX or Pure Charge&Go BCT IX
2) Not available in Insio ITC/ITE IX
3) Not available in Motion Charge&Go SP IX
†) Motion Charge&Go IX BTE family only available in performance level 7IX, 5IX, 3IX
††) Two-microphone models only

Features	Description	Performance levels: Styletto IX, Pure Charge&Go IX family, Motion Charge&Go IX family†, Insio IX ITC††/ITE IX					Performance levels: Active Pro IX	
		7 IX	5 IX	3 IX	2 IX	1 IX	Active Pro	Active
Fitting Optimization								
Processing Channels / Gain & MPO handles	Number of processing channels/handles in Connexx	48/20[1]	36/18	32/16	24/12	16/8	48/20	16/8
HD Music	Preset optimized music programs for enhanced enjoyment of non-streamed music.	3	3	1	1		3	
Hearing programs	Number of programs configurable in Connexx.	6	6	6	4	4	6	4
CROS/BiCROS compatibility	A CROS device worn on an ear with unaidable hearing loss transmits sound to a hearing aid worn on the better ear.	✓	✓	✓	✓	✓	✓	✓
Frequency Compression	Improves the detection of speech sounds for wearers with reduced high frequency hearing thresholds. Frequency compression is selectively activated based on audiological criteria.	✓	✓	✓	✓	✓	✓	✓
Smart Optimizer and Data Logging	Tailored fine-tuning recommendations based on logging and analysis of the wearer's hearing aid use.	✓	✓	✓	✓	✓	✓	✓
Acclimatization Manager	Allows the wearer to benefit from enhanced audibility in the longer term, while enjoying optimal device acceptance when initially fitted.	✓	✓	✓	✓	✓	✓	✓
InSituGram	Allows quick and convenient measurement of hearing thresholds with the hearing aids in place.	✓	✓	✓	✓	✓	✓	✓
AutoFit™	Automatically provides an optimized real-ear insertion gain (REIG) setting without the need to manually match targets, enabling significantly faster real-ear verification.	✓	✓	✓	✓	✓	✓	✓

1) Motion Charge&Go SP IX has 36/18 channels/handles
†) Motion Charge&Go IX BTE family only availble in performance level 7IX, 5IX, 3IX
††) Two-microphone models only

Feature overview

Silk Charge&Go IX, Insio Charge&Go CIC IX, Insio IIC/CIC/ITC† IX

Features	Description	Performance levels				
		7 IX	5 IX	3 IX	2 IX	1 IX
Conversation management						
Dynamic Soundscape Processing 3.0	Integrates inputs for precise hearing aid optimization, including (when available): acoustic sensors, noise detection and SNR analysis.	●●●●●	●●●●○	●●●○○	●●○○○	●●○○○
Binaural OneMic Directionality 2.0	Microphones from each hearing aid work together to create a directional listening pattern preferentially amplifying speech from the front to support hearing in noisy environments.[1]	✓	✓	✓		
SpeechStabilizer	Preserves the natural envelope of speech in any situation, including very noisy situations, which preserves acoustic cues that are important for speech understanding.	✓	✓	✓	✓	✓
Speech and Noise Management	Reduces listening effort using sophisticated sound processing.	✓	✓	✓	✓	✓
e2e wireless 4.0	Controls coupling, synchronization, CROS/BiCROS, and a range of sound features.[1]	✓	✓	✓	✓	✓
Listening Comfort						
Extended bandwidth	12 kHz available bandwidth for enhanced processing of high-frequency speech and environmental sounds.	✓				
Auto EchoShield	Automatically reduces the effect of reverberations and sudden loud noises in reverberant environments.	✓				
Reverberant Room Program	Dedicated program for reverberant environments.	✓	✓			
SoundSmoothing™	Reduces the annoyance of sudden loud impact sound such as when handling glassware.	●●●●●	●●●●●	●●●●●	●●●○○	●○○○○
eWindScreen™	Reduces the interference of wind noise when outdoors.	●●●●●	●●●●●	●●●○○	●○○○○	●○○○○
Extended dynamic range	117 dB 24 bit processing for enhanced digital encoding of soft and loud sounds.	✓	✓	✓	✓	✓
Feedback Cancellation	High speed monitoring and control of feedback in individual processing channels.	✓	✓	✓	✓	✓
Wearer Empowerment & Connectivity						
Signia App (iOS and Android)	Easy-to-use smartphone app enables remote control adjustments and TeleCare.	✓	✓	✓	✓	✓
TeleCare	Provides remote services including the TeleCare Dashboard and CareChat (convenient text messages, phone calls, and video calls).	✓	✓	✓	✓	✓
miniPocket compatibility	Optional keychain remote control.	✓	✓	✓	✓	✓
Telecoil	Available as an option for Insio ITC IX. Enables direct connection to audio loop systems and compatible telephone signals.	✓	✓	✓	✓	✓

1) Available for bilateral fittings with e2e wireless link
†) One-microphone models

Features	Description	Performance levels				
		7 IX	5 IX	3 IX	2 IX	1 IX
Tinnitus						
Ocean Waves Therapy Signal (presets)	Therapy signals that mimic the sound of the sea to provide a soothing listening experience and reduce the impact of tinnitus.	4	4	4	4	
Static Therapy Signal (channels / presets)	Soft, gentle therapy signals that reduce the impact of tinnitus.	20/5	18/5	16/5	12/5	
Notch Therapy	Signia's unique, patented tinnitus therapy which is proven especially effective against tonal tinnitus.	✓	✓	✓	✓	
Fitting Optimization						
Processing Channels / Gain & MPO handles	Number of processing channels/handles in Connexx	48/20	36/18	32/16	24/12	16/8
HD Music	Preset optimized music programs for enhanced enjoyment of non-streamed music.	3	3	1	1	
Hearing programs	Number of programs configurable in Connexx.	6	6	6	4	4
CROS/BiCROS compatibility	A CROS device worn on an ear with unaidable hearing loss transmits sound to a hearing aid worn on the better ear.[2]	✓	✓	✓	✓	✓
Frequency Compression	Improves the detection of speech sounds for wearers with reduced high frequency hearing thresholds. Frequency compression is selectively activated based on audiological criteria.	✓	✓	✓	✓	✓
Smart Optimizer and Data Logging	Tailored fine-tuning recommendations based on logging and analysis of the wearer's hearing aid use.	✓	✓	✓	✓	✓
Acclimatization Manager	Allows the wearer to benefit from enhanced audibility in the longer term, while enjoying optimal device acceptance when initially fitted.	✓	✓	✓	✓	✓
InSituGram	Allows quick and convenient measurement of hearing thresholds with the hearing aids in place.	✓	✓	✓	✓	✓
AutoFit™	Automatically provides an optimized real-ear insertion gain (REIG) setting without the need to manually match targets, enabling significantly faster real-ear verification.	✓	✓	✓	✓	✓

1) Available for bilateral fittings with e2e wireless link
2) Requires a hearing aid with e2e wireless link

Signia AX

Revolutionary split processing

Performance levels

- 7AX
- 5AX
- 3AX
- 2AX
- 1AX

Hearing isn't always easy. A big group of people talking at the same time, softly spoken talkers, too much background noise… Sometimes sounds blend together making it impossible to focus on what you want to hear.

AX Augmented Xperience

Augmented Xperience changes the way your clients hear the world. Instead of them having to strain to discern speech from background noise, our new platform uses proprietary Augmented Focus™ technology to process the two separately and create a clear contrast. It then recombines them to deliver outstanding speech clarity in a fully immersive environment.
For instance, 100% of participants in a recent study[1] reported excellent speech understanding in their home environment using this technology.

Augmented Xperience is available across all performance levels from 7 to 1 so you can offer all your customers outstanding speech clarity in an immersive environment

ITE



Insio Charge&Go AX

RIC



Pure Charge&Go AX



Pure Charge&Go T AX



Pure 312 AX

SLIM-RIC



Styletto AX

(1) Jensen, Pischel, Taylor, Schulte: Performance of Signia AX in At-Home Listening Situations (Signia White Paper, 2021)



Designed for brilliant hearing and looks



Key features

- Augmented Focus
- Own Voice Processing 2.0
- Signia Assistant
- Android & iPhone connectivity
- Recharge on-the-go
- CROS compatibility
- Award-winning design*

Performance levels

7AX
5AX
3AX
2AX
1AX

sDemo

Styletto AX

Our new Styletto AX rechargeable hearing aid looks as amazing as its sounds. Its award-winning SLIM-RIC design* is smaller than before for an even more discreet fit. Built on our leading-edge Augmented Xperience platform, Styletto AX uses our proprietary Augmented Focus™ technology to split speech from surrounding sounds, process them separately, and recombine them to deliver outstanding speech clarity in an immersive environment.

Also available as a CROS Styletto AX.



Charging options
Styletto AX offers the choice of a portable charger with Qi technology – wearers simply place the charger on a charging pad instead of having to plug it into the mains. Alternatively, they can opt for a Dry&Clean charger, which simultaneously charges the hearing aids and dries and cleans them via UV light. See pages 60-61 for more information on chargers.

SLIM-RIC

Model			
Receiver type	miniReceiver 3.0 S	miniReceiver 3.0 M	miniReceiver 3.0 P
			
Performance			
Matrix (2 ccm)	46/110 dB SPL	60/119 dB SPL	65/122 dB SPL
Battery			
Type	Lithium-ion power cell	Lithium-ion power cell	Lithium-ion power cell

* iF Design Award 2019 for Styletto

Immersive sound in our smallest rechargeable RIC



Key features

- Augmented Focus
- Own Voice Processing 2.0
- Signia Assistant
- Android & iPhone connectivity
- Convenient rechargeability
- CROS compatibility
- Available as T-Coil Option

Performance levels

7AX
5AX
3AX
2AX
1AX

sDemo

Pure Charge&Go AX / Pure Charge&Go T AX

Now available in all performance levels, Pure Charge&Go AX is ideal for experienced and first-time hearing aid wearers who are looking for a sleek aesthetic device that is fully featured. It combines an attractive design with proven rechargeability, state-of-the-art Bluetooth connectivity for Android and iPhone, and the landmark AI of the Signia Assistant.

Also available as a CROS Pure Charge&Go AX, or as Pure Charge&Go T AX with integrated T-Coil.



Charging options
The Pure Portable Charger offers comfort and charging on the go. The Pure Dry&Clean Charger also offers UV cleaning and drying to keep the hearing aids in peak condition. And the small, reliable standard Pure Charger meets basic desktop charging needs. Also compatible with the Multi Charger. See pages 60-61 for more information on chargers.

RIC

Model				
Receiver type	miniReceiver 3.0 S	miniReceiver 3.0 M	miniReceiver 3.0 P	miniReceiver 3.0 HP
				
Performance				
Matrix (2 ccm)	46/110 dB SPL	60/119 dB SPL	65/122 dB SPL	75/131 dB SPL
Battery				
Type	Lithium-ion power cell	Lithium-ion power cell	Lithium-ion power cell	Lithium-ion power cell

The first custom hearing aids with contactless charging and Bluetooth

Key features

- Augmented Focus
- Signia Assistant
- Android & iPhone connectivity
- Contactless charging
- Custom fit
- Award-winning design*



Performance levels

7AX

5AX

3AX

2AX

1AX

Insio Charge&Go AX

Our new rechargeable Insio Charge&Go AX ITC and ITE models redefine what custom hearing aids can offer experienced wearers. The smaller, more discreet design offers a higher fitting rate than previous devices and all-day comfort, while our leading-edge Augmented Xperience platform delivers outstanding speech clarity in an immersive environment. The custom hearing aids even include a motion sensor to adapt their sound optimally to the wearer's situation by detecting if they are moving or stationary.



Insio Charger
These are the first custom hearing aids with contactless charging. Wearers simply place Insio Charge&Go AX in the Insio Charger without having to worry about any exact charging contact.

Model	ITC/HS ITE	ITC/HS ITE	ITC/HS
			
Performance			
Matrix (2 ccm)	124/65 dB SPL	118/55 dB SPL	113/50 dB SPL
Battery			
Type	Lithium-ion power cell		

* CES Innovation awards 2022 honoree and iF Design Award 2022 for Insio Charge&Go AX

A RIC for outstanding speech clarity

Key features

- Augmented Focus
- Own Voice Processing 2.0
- Signia Assistant
- Android & iPhone connectivity
- CROS compatibility



Performance levels

7AX

5AX

3AX

2AX

1AX

sDemo

Pure 312 AX

Equipped with all the benefits of our Augmented Xperience platform for outstanding speech clarity, Pure 312 AX and CROS Pure 312 AX attract wearers with their fresh, elegantly slim design. Despite their discreet form, these tiny hearing aids provide fully-featured Bluetooth connectivity for maximum enjoyment of music, calls and TV from Android and iOS devices, plus all the benefits of our Signia Assistant. Pure 312 AX's Zn-Air battery delivers a runtime of around 4 days.*

Model	RIC			
Receiver type	miniReceiver 3.0 S	miniReceiver 3.0 M	miniReceiver 3.0 P	miniReceiver 3.0 HP
				
Performance				
Matrix (2 ccm)	44/108 dB SPL	58/116 dB SPL	62/119 dB SPL	74/128 dB SPL
Battery				
Size	312	312	312	312
Power consumption	1.5mA	1.7mA	1.7mA	1.7mA

*Based on 17 hours wearing per day including 5 hours of streaming

Augmented Xperience feature overview

Our leading-edge Augmented Xperience platform is the first to offer outstanding speech clarity in an immersive environment, combined with connectivity for Android[1] and iOS devices plus artificial intelligence (AI) for live remote assistance.

Features	Information	Performance levels				
		7AX	5AX	3AX	2AX	1AX
Improved AX Soundscape Processing	Dynamic Soundscape Processing takes sound processing to a more precise level by including motion data for a more realistic sound adapted to everyday life.	●●●●●	●●●●●	●●●●●	●●●○○	●●●○○
Augmented Focus™	Separates the focus sound (speech) from surroundings, creating a clear contrast between the two, which makes it easier to follow conversations.	✓	✓	✓	✓	✓
Acoustic Sensors	The redesigned acoustic sensors can detect more variables of sound than ever before, providing more accurate information about the environment.	✓	✓	✓	✓	✓
e2e wireless 4.0	Control coupling, synchronization, binaural directionality, Spatial Speech Focus, CROS/ BiCROS.	✓	✓	✓	✓	✓
Motion Sensor	An integrated motion sensor detects whether the wearer is walking, running or sitting still. This motion data complements the acoustic data to optimize Dynamic Soundscape Processing 2.0.	✓	✓	✓		
OVP2.0 [3]	OVP2.0 utilizes real-time recognition of the wearer's voice to deliver a natural own-voice impression. Available for RIC/BTE devices with e2e wireless in binaural fittings.	✓	✓	✓		
Sound Clarity	Sound Clarity provides lifelike sound which covers the full spectrum of wearing situations.	●●●●●	●●●●○	●●●○○	●●●○○	●●○○○
Signal processing / Gain & MPO	Channels/handles in Connexx. For Signia AX, channels are per processing chain (Focus/Surroundings).	48/20	32/16	24/12	16/8	16/8
Hearing Programs	Number of programs configurable in Connexx.	6	6	6	4	4
Extended dynamic range	117 dB 24 bit processing for enhanced digital encoding of soft and loud sounds.	✓	✓	✓	✓	✓
Speech and noise management	Reduces listening effort using sophisticated sound processing.	✓	✓	✓	✓	✓
SoundSmoothing™	Reduces the annoyance of sudden loud impact sound such as when handling glassware.	✓	✓	✓	✓	
Feedback cancellation	High speed monitoring and control of feedback in individual processing channels.	✓	✓	✓	✓	✓
HD Music	Preset programs for enhanced enjoyment of non-streamed music.	3	3	1	1	
eWindScreen	Reduces the annoyance of wind noise when outdoors.	✓	✓	✓	✓	
Extended bandwidth	12 kHz bandwidth for enhanced processing of high-frequency speech and environmental sounds.	✓				
EchoShield	Dedicated program for reverberant environments.	✓	✓			
Speech Quality	Speech Quality delivers improved speech understanding and communication ease in noise.	●●●●●	●●●●○	●●●○○	●●●○○	●●●○○
Binaural Directionality	Directionality provides scientifically proven enhanced speech understanding in noise. Available for bilateral fittings with directional microphones and e2e wireless link.	✓	✓	✓		
Wireless CROS/BiCROS	Requires CROS transmitters.	✓	✓	✓	✓	✓
Frequency compression	Improves the detection of speech sounds for wearers with reduced high frequency hearing thresholds. Frequency compression is selectively activated based on audiological criteria.	✓	✓	✓	✓	✓

Features	Information	Performance levels				
		7AX	5AX	3AX	2AX	1AX
Spatial SpeechFocus	Depending on the direction of the dominant speech source, Spatial SpeechFocus steers the directional beam to the left, right, or behind the wearer automatically in a car situation, or in a dedicated "Strolling" program. Available for bilateral fittings with directional microphones and e2e wireless link.	✓1)	✓2)			
Wearer Interaction	Technically inclined wearers may appreciate having extra control over their hearing aids. Signia Augmented Xperience hearing aids offer these wearers a comprehensive range of adjustment options via the Signia app.	●●●●●	●●●●●	●●●●●	●●●●●	●●●●●
Signia Assistant	Signia Assistant offers wearers instant 24/7 support and troubleshooting around the clock. It provides detailed information in Connexx to address each client's fine-tuning needs more precisely.	✓	✓	✓	✓	✓
Signia app (iOS and Android)	The new Signia app combines all functionalities of previous apps (touchControl, myControl, myHearing app). The functionality offered in the app is automatically customized depending on the connected hearing aid.	✓	✓	✓	✓	✓
Adaptive Streaming Volume	Speech intelligibility weighted adaptation of the streaming volume depending on ambient noise.	✓	✓	✓	✓	✓
Spatial Configurator	Spatial Configurator is part of the Signia app and allows the wearer to adjust the direction of the microphone beam.	✓	✓			
Direct streaming for iOS	Signia Augmented Xperience hearing aids can be connected directly to Apple devices for phone calls and audio streaming – without additional accessories.	✓	✓	✓	✓	✓
Direct streaming for Android	Direct audio streaming from ASHA compatible Android devices to Bluetooth enabled Augmented Xperience hearing aids.	✓	✓	✓	✓	✓
StreamLine TV accessory	Stream high quality stereo sound directly from a TV to Bluetooth equipped Augmented Xperience hearing aids.	✓	✓	✓	✓	✓
Call Control/ Tap Control + HandsFree for iOS	Call Control/Tap Control: accept and end call using a hearing aid Handsfree for iOS: hands-free mobile phone calls – even in noisy environments, combine phone calls with other activities or to listen to music or podcasts via your hearing aids	✓	✓	✓		
StreamLine Mic accessory	Stream high quality stereo sound directly from any Bluetooth device (Android devices, older Apple devices, other smartphones, PCs, laptops, etc.) to Bluetooth equipped Augmented Xperience hearing aids.	✓	✓	✓	✓	✓
Tinnitus	A choice of tinnitus treatments based on amplified notch therapy or traditional noise therapy signals. It is even possible to combine the two approaches.	✓	✓	✓	✓	
Notched Amplification Therapy	Innovative tinnitus therapy delivered via everyday amplification.	✓	✓	✓	✓	
Smart Optimizer and Data Logging	Smart Optimizer (found under Data Logging) provides tailored fine-tuning recommendations based on an analysis of the wearer's usage of the hearing aid.	✓	✓	✓	✓	✓
Acclimatization manager	Allows the wearer to benefit from enhanced audibility in the longer term, while enjoying optimal device acceptance when initially fitted.	✓	✓	✓	✓	✓
InSituGram	Allows quick and convenient screening audiograms and checking of thresholds while the hearing aids are worn.	✓	✓	✓	✓	✓
AutoFit	Automatically provides an optimized real-ear insertion gain (REIG) setting without the need to manually match targets. The new IMC2 standard expands the range of REM equipment which supports AutoFit with Connexx and Noah 4.5.1 or higher.	✓	✓	✓	✓	✓
TeleCare	TeleCare provides remote services including additional tools to follow-up with patients. In addition, user engagement and autonomy can be increased via the Signia app.	✓	✓	✓	✓	✓
Remote Services	Includes TeleCare Dashboard, Remote Fine Tuning, CareChat (incl. text messages, phone calls and video calls), virtual follow-up.	✓	✓	✓	✓	✓

1) In Universal program: Automatic mic beam to the right / left / rear directions.
2) In Universal program: Automatic mic beam to the rear, surround mode for side speech. In manual programs, microphone beam to the right/left/rear directions.
3) Not available for Insio Charge&Go AX

Signia Xperience

Empowering people to hear what matters

Performance levels

7X
5X
3X
2X
1X

Features	Information	Performance levels				
		7X	5X	3X	2X	1X
YourSound Technology	Redesigned acoustic sensors plus an integrated motion sensor provide detailed information about where the wearers are and what they are doing. Dynamic Soundscape Processing and Own Voice Processing (OVP™) process this information delivering the most clear and personalized sound.	●●●●●	●●●●●	●●●●●	●●●○○	●●●○○
Acoustic Sensors	Detects more variables of sound than ever before, providing more accurate information about the environment.	✓	✓	✓	✓	✓
Dynamic Soundscape Processing	Steers Sound Clarity and Sound Quality features for natural sound and speech in every situation from any direction, even when moving.	●●●●●	●●●●●	●●●●●	●●●○○	●●●○○
Motion Sensor[A]	Detects whether the wearer is walking, running or not. This information is used for Dynamic Soundscape Processing.	✓	✓	✓		
OVP™ (Own Voice Processing)[A]	Utilizes real-time recognition of the wearer's voice to deliver a natural own-voice impression.	●●●●●	●●●●●	●●●●●		
Signal Processing / Gain & MPO	Channels/handles in Connexx.	48/20[B]	32/16	24/12	16/8	16/8
Signia Assistant[A]	Offers the wearers instant 24/7 support and troubleshooting around the clock. It provides detailed information in Connexx to address each client's fine-tuning needs more precisely.	✓	✓	✓	✓	✓
TeleCare	Provides virtual appointments, remote fine-tuning, progress updates, hearing lessons and daily satisfaction ratings.	✓	✓	✓	✓	✓
Tinnitus	A choice of tinnitus treatments based on amplified notch therapy or traditional noise therapy signals. It is even possible to combine the two approaches.	✓	✓	✓	✓	
eWindScreen™	Reduces the annoyance of wind noise when outdoors, improving sound comfort. eWindScreen Binaural requires an e2e wireless link.	Binaural[A] auto	Binaural[A] manual	Monaural	Monaural	
Binaural OneMic Directionality	Delivers a highly sophisticated speech-in-noise performance. Available for Silk X.	✓	✓	✓		
Spatial SpeechFocus	Detects the loudest voice in the environment and adjusts the directional beam towards that voice (i.e. front, side, or backward facing). Available for bilateral fittings with directional microphones and e2e wireless link.	front,back, right,left	front,back, omni			
Extended Bandwidth	12 kHz bandwidth for enhanced processing of high-frequency speech and environmental sounds.	✓				
EchoShield	Dedicated program for reverberant environments.	✓				

A) Not available with Silk X
B) For Motion Charge&Go SP X: 34/16

Experience the full potential of BTEs

Performance levels

7X
5X
3X
2X
1X

sDemo



Key features

- Most natural own voice and best speech understanding in noise, even when the wearer is in motion
- World's first integrated acoustic-motion sensors
- Long-lasting high quality audio streaming
- TeleCare
- Remote Control
- Connectivity & Streaming with StreamLine Mic and StreamLine TV

Signia app ready

Smaller Design

Equipped with Bluetooth

Rechargeable

OVP Own Voice Processing



Inductive Charger II
This charger comes with a charger lid to protect the hearing aids, a fast charging option, and a drying function. Its smart design accommodates custom molds.

Motion Charge&Go X

Comparable in its small size to a 312- device, Motion Charge&Go X is the ultimate BTE combining sleek design with rechargeability, Bluetooth, and the Signia Assistant. It is flexible to fit, easy to use and comfortable to wear, making this an easy choice for those looking for convenience and robustness – all in one.

BTE

Model	Earhook	ThinTube 3.0	ThinTube 3.0 P
			
Performance			
Matrix (2 ccm)	60/133 dB SPL	60/125 dB SPL	63/126 dB SPL
Battery			
Type	Lithium-ion power cell	Lithium-ion power cell	Lithium-ion power cell

Hear what really matters

Performance levels

7X

5X

3X

2X

1X

sDemo



Motion Charge&Go P X

Smaller than its predecessor without compromising on power, Motion Charge&Go P X combines Signia Xperience hearing technology with the complete convenience of Li-ion recharging, making annoying battery changes a thing of the past. This powerful hearing solution for every situation with T-Coil ensures that whatever you spend your day doing, you will be able to perform at your best.

Key features

- Most natural own voice and best speech understanding in noise, even when the wearer is in motion
- World's first integrated acoustic-motion sensors
- Long-lasting high quality audio streaming
- TeleCare
- Remote Control
- Connectivity & Streaming with StreamLine Mic and StreamLine TV

Signia app ready

Smaller Design

Equipped with Bluetooth

Rechargeable

Integrated T-Coil

OVP Own Voice Processing



Charging options
The Motion P X Dry&Clean Charger also offers UV cleaning and drying to keep the hearing aids in peak condition. And the small, reliable standard Motion P X Charger meets basic desktop charging needs. Motion Charge&Go P X is also compatible with the Multi Charger.

BTE

Model	Earhook	ThinTube 3.0	ThinTube 3.0 P
			
Performance			
Matrix (2 ccm)	77/135 dB SPL	66/130 dB SPL	70/131 dB SPL
Battery			
Type	Lithium-ion power cell	Lithium-ion power cell	Lithium-ion power cell

Superpower your hearing

Performance levels

7X

5X

3X

2X

1X

sDemo



Motion Charge&Go SP X

This is the world's first rechargeable Super Power BTE. Closing the gap between Super Power and super features, the innovative Motion Charge&Go SP X ultra-power device is the complete package for those who want to Be Brilliant despite severe to profound hearing loss. With Bluetooth connectivity, Li-ion rechargeability and T-Coil, it offers a world of new hearing and handling possibilities not familiar to those with conventional SP devices.

Key features

- Most natural own voice and best speech understanding in noise, even when the wearer is in motion
- World's first integrated acoustic-motion sensors
- Long-lasting high quality audio streaming
- TeleCare
- Remote Control
- Connectivity & Streaming with StreamLine Mic and StreamLine TV

Signia app ready

Smaller Design

Equipped with Bluetooth

Rechargeable

Integrated T-Coil

OVP Own Voice Processing



Charging options
The Motion SP X Dry&Clean Charger also offers UV cleaning and drying to keep the hearing aids in peak condition. And the small, reliable standard Motion SP X Charger meets basic desktop charging needs. Motion Charge&Go SP X is also compatible with the Multi Charger.

BTE

Model	Earhook
	
Performance	
Matrix (2 ccm)	82/140 dB SPL
Battery	
Type	Lithium-ion power cell

Discretion guaranteed

Performance levels

7X
5X
3X
2X
1X

sDemo



Silk X

Silk X is the ideal solution for those who value discretion above all. The only ready-to-wear CIC available on the Signia Xperience platform delivers clear speech understanding in an ultra-discreet size with instant fit. It is tailor-made without the wait: Your clients can walk in, experience Silk X and walk out the door with a new CIC without having to wait for a custom-built device. Thanks to its flexible silicone Click Sleeves, Silk X sits securely in the ear for a brilliant fit on the spot.

Key features


- Nearly invisible in the ear for complete discretion
- Signia Xperience technology for a superior hearing experience
- Click Sleeves for a quick and easy instant-fit
- Remote control
- TeleCare

Signia app ready

Compatible with CROS device for people with unaidable hearing loss on one side

Instant fit

DSP Dynamic Soundscape Processing

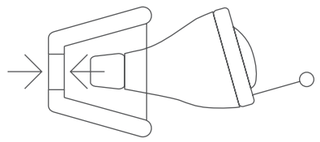

One click and you're done.
The sleeves come in four different sizes and a closed and open build version suitable for most ear anatomies and hearing needs.

Instant Fit ITE

Model	Click CIC
	
Performance	
Matrix (2 ccm)	50/114 dB SPL
Battery	
Size	10
Power consumption	1.2 mA

Fully versatile, fully connected

Performance levels

7X
5X
3X
2X
1X

sDemo



Pure 312 X

Equipped with all the benefits of the Signia Xperience platform, Pure 312 X attracts wearers with its elegantly slim design. Despite its discreet form, Pure 312 X provides fully-featured Bluetooth connectivity for maximum enjoyment of music, calls and TV.

An integrated T-Coil is also available.

Key features

- Most natural own voice and best speech understanding in noise, even when the wearer is in motion
- World's first integrated acoustic-motion sensors
- Long-lasting high quality audio streaming
- TeleCare
- Remote Control
- Connectivity & Streaming

Signia app ready

Equipped with Bluetooth

Optional T-Coil

Compatible with CROS device for people with unaidable hearing loss on one side

OVP Own Voice Processing

RIC

Model				
Receiver type	miniReceiver 2.0 S	miniReceiver 2.0 M	miniReceiver 2.0 P	miniReceiver 2.0 HP
				
Performance				
Matrix (2 ccm)	45/108 dB SPL	60/119dB SPL	70/124 dB SPL	75/130 dB SPL
Battery				
Size	312	312	312	312
Power consumption	1.7 mA	1.9 mA	1.8 mA	1.8 mA

Sirion Connect

Connect with your expectations

With Sirion Connect, we're once again redefining what's possible in the essentials category, bringing seamless connectivity, exceptional design and advanced technology to more people than ever before.

Effortless to use with direct Bluetooth® connectivity via MFi or ASHA, tailored control via the Signia app, and a streamlined wireless fitting calculated in 48 channels for first fit precision.

Elevated design with convenient tactile push buttons, a premium finish, and a robust IP68 rating.

Enhanced technology with Signia's unique split processing technology, which processes speech and noise signals separately for prominent speech within a realistic soundscape.

RIC



Sirion Connect RIC

BTE



Sirion Connect P



Sirion Connect SP

Connect with your expectations



Key features

- Augmented Focus
- e2e Wireless 4.0
- Bluetooth connectivity
- TeleCare
- Signia app

Sirion Connect

Bringing connectivity and Signia's unique split processing technology, Augmented Focus, to more people than ever.

	RIC			
Model				
Receiver type	S Receiver	M Receiver	P Receiver	HP Receiver
				
Performance				
Matrix (2 ccm)	46/109 dB SPL	58/117 dB SPL	63/120 dB SPL	74/128 dB SPL
Battery (Zinc-Air)				
Size	312	312	312	312

Connect with your expectations



Sirion Connect P



Sirion Connect SP

Sirion Connect P

Sirion Connect SP

Bringing connectivity and Signia's unique split processing technology, Augmented Focus, to more people than ever. Two BTE power levels - including a superpower device - cover the full range of amplification needs.

Key features

- Augmented Focus
- e2e Wireless 4.0
- Bluetooth connectivity
- TeleCare
- Signia app
- Integrated telecoil

	BTE P			BTE SP
Model				
	Earhook	ThinTube 3.0	ThinTube P 3.0	Earhook
				
Performance				
Matrix (2 ccm)	77/135 dB SPL	66/129 dB SPL	70/130 dB SPL	82/143 dB SPL
Battery (Zinc-Air)				
Size	13	13	13	675

Sirion Connect feature overview

Features	Model		
	RIC	BTE P	BTE SP
Processing Channels / Gain & MPO handles	12/6	12/6	12/6
Hearing programs	4	4	4
Augmented focus	●	●	●
Speech and Noise Management	●	●	●
Feedback cancellation	●	●	●
e2e Wireless 4.0	●	●	●
TruEar 360	●	●	●
Automatic directionality	●	●	●
InSituGram	●	●	●
AutoFit™	●	●	●
Direct streaming for iOS (MFi) & Android (ASHA)	●	●	●
Signia App (iOS and Android)	●	●	●
Streamline TV compatibility	●	●	●
StreamLine Mic compatibility	●	●	●
miniPocket compatibility	●	●	●
Telecoil	-	●	●

Orion Charge&Go

Effortless to use
Elevated design
Enhanced technology

Wearers in the essentials category can say goodbye to the hassle of small batteries and hello to the convenience of rechargeability.

With a premium finish, Orion Charge&Go hearing aids are not only good-looking but also hard-wearing, with an IP68 rating for resistance to dust and moisture. And with intelligent hearing technology that adapts to the sound environment, wearers can enjoy sophisticated hearing support wherever they are.

Benefit from an easy wireless fitting process, with a first fit calculated in 48 channels ensuring an efficient, personalized fit right from the start.

RIC



Orion Charge&Go RIC

BTE



Orion Charge&Go P
Orion Charge&Go P T



Orion Charge&Go SP
Orion Charge&Go SP T

Charge your expectations

Performance levels

- 200
- 100
- 75
- 50
- sDemo



Orion Charge&Go RIC

An attractive and comfortable RIC with four receiver options to cover a wide range of amplification needs.

Key features

- TruEar
- Adaptive directionality
- eWindScreen
- Convenient rechargeability
- e2e Wireless
- Signia app

	RIC			
Model				
Receiver type	miniReceiver 2.0 S	miniReceiver 2.0 M	miniReceiver 2.0 P	miniReceiver 2.0 HP
				
Performance				
Matrix (2 ccm)	45/108 dB SPL	60/119 dB SPL	70/124 dB SPL	75/130 dB SPL
Battery				
Type	Lithium-ion power cell	Lithium-ion power cell	Lithium-ion power cell	Lithium-ion power cell

Charge your expectations



Orion Charge&Go P
Orion Charge&Go P T

Orion Charge&Go SP
Orion Charge&Go SP T

Key features

- TruEar
- Speech and noise management
- eWindScreen
- Convenient rechargeability
- e2e Wireless
- Signia app
- Optional telecoil

Orion Charge&Go P

Performance levels

200
100
75
50

Orion Charge&Go SP

Orion Charge&Go P T

Performance levels

200
100
sDemo

Orion Charge&Go SP T

Orion Charge&Go BTEs are available in two power levels - P and SP - to cover the full range of amplification needs. Orion Charge&Go P T and Orion Charge&Go SP T offer integrated telecoil.

BTEs

Model	Orion Charge&Go P Orion Charge&Go P T	Orion Charge&Go SP Orion Charge&Go SP T
		
Performance		
Matrix (2 ccm)	70/134 dB SPL	80/140 dB SPL
Battery		
Type	Lithium-ion power cell	Lithium-ion power cell

Orion Charge&Go feature overview

Features	Performance levels			
	200	100	75	50
Frequency channels / compression controls	16/8	16/8	12/6	8/4
Number of programs	4	4	4	4
48-channel first fit	●	●	●	●
Adaptive compression	●	●	●	●
Extended dynamic range	●	●	●	●
Signia app	●	●	●	●
miniPocket compatibility	●	●	●	●
InSituGram	●	●	●	●
AutoFit	●	●	●	●
Smart optimizer and Data logging	●	●	●	●
Feedback cancellation	●	●	●	●
Speech and noise management	●	●	●	●
Microphone directionality	●	●	●	●
Automatic directionality	●	●	●	-
Adaptive directionality	●	-	-	-
TruEar	●	●	-	-
SoundBalance	●	●	-	-
Frequency compression	●	●	-	-
Automatic acclimatization manager	●	●	-	-
e2e wireless	●	●	-	-
User control coupling via e2e	●	●	-	-
SoundSmoothing™	●	-	-	-
eWindScreen™	●	-	-	-
Tinnitus static therapy signal	●	-	-	-

Intuis 4

Advanced technology for more of your clients

Advanced technology for more of your clients

Key features

- Augmented Focus
- Own Voice Processing 2.0
- Equipped with Bluetooth
- TeleCare
- Signia app
- Integrated telecoil



Intuis 4 BTEs

Intuis 4 BTE hearing aids feature advanced technology that delivers brilliant sound and speech clarity, plus the convenience of enhanced connectivity and Bluetooth streaming. Suitable for all levels of hearing loss, the ergonomic design and easy-to-use controls of Intuis 4 make daily operation straightforward and comfortable. The Signia app provides a suite of controls and communication options.

Features	Intuis 4					
	4.7	4.5	4.3	4.2	4.1	4.0
Processing Channels / Gain & MPO handles	48/20[1]	32/16	24/12	16/8	16/8	16/6
Hearing programs	6	6	6	4	4	4
Augmented Focus	●	●	●	●	●	●
Speech and Noise Management	●	●	●	●	●	●
Feedback Cancellation	●	●	●	●	●	●
e2e wireless 4.0	●	●	●	●	●	●
InSituGram	●	●	●	●	●	●
AutoFit™	●	●	●	●	●	●
miniPocket compatibility	●	●	●	●	●	●
Signia App (iOS and Android)	●	●	●	●	●	●[2]
Signia Assistant	●	●	●	●	●	-
Direct streaming for iOS (MFi) & Android (ASHA)	●	●	●	●	●	-
Streamline TV and StreamLine Mic	●	●	●	●	●	-
TruEar 360	●	●	●	●	●	-
Automatic Directionality	●	●	●	●	●	●
Adaptive Directionality	●	●	●	●	●	-
Binaural Directionality	●	●	●	-	-	-
SpeechFocus	●	●	-	-	-	-
Spatial SpeechFocus	●	-	-	-	-	-
CROS/BiCROS compatibility	●	●	●	●	●	-
Frequency Compression	●	●	●	●	●	-
Acclimatization Manager	●	●	●	●	●	-
Tinnitus therapy	●	●	●	●	-	-
SoundSmoothing™	IIIII	IIIII	IIIII	III	-	-
eWindScreen™	IIIII	IIIII	III	I	-	-
Own Voice Processing 2.0 (OVP™ 2.0)	●	●	●	-	-	-
Spatial Configurator	●	●	-	-	-	-
Reverberant Room Program	●	●	-	-	-	-
Auto EchoShield	●	-	-	-	-	-
Extended bandwidth	●[3]	-	-	-	-	-

[1] 34/16 for Intuis SP 4
[2] Acoustic pairing only, Signia App with reduced functionality.
[3] Intuis M4 and Intuis P 4 only

	BTE		
Model			
	Intuis M 4	Intuis P 4	Intuis SP 4
			
Performance			
Matrix (2 ccm)	ThinTube 3.0: 61/125 dB SPL ThinTube 3.0 P: 64/126 dB SPL Earhook: 60/133 dB SPL	ThinTube 3.0: 66/129 dB SPL ThinTube 3.0 P: 70/130 dB SPL Earhook: 77/135 dB SPL	Earhook: 82/140 dB SPL
Battery			
Size	13	13	675

Intuis 3

Proven technology, outstanding comfort

Whether it's a receiver-in-canal (RIC), a behind the ear (BTE), or an in-the-ear (ITE) model, the Intuis 3 family has a solution for you. Suitable for mild to severe hearing loss, all hearing aids provide you with a reliable, clear and comfortable listening experience and a maximum of discretion with remote control options via the Signia app or the miniPocket™.

Feature overview

Features	Intuis 3
Frequency channels / compression controls	12/6
Number of programs	4
Speech and noise management	●
Feedback cancellation	●
Directional microphones	automatic[B]
e2e wireless	●
Signia app	●
Smart optimizer and Data logging	●
InSituGram	●
AutoFit	●
miniPocket	●

B) Only BTE

RIC



Intuis 3 RIC312

BTE



Intuis 3 S



Intuis 3 M



Intuis 3 P



Intuis 3 SP

ITE



Intuis 3 Click CIC, Click ITC



Intuis 3 CIC, ITC, ITE

Small and powerful



Key features

- For mild to severe hearing loss
- Small design for a high level of wearing comfort
- TeleCare
- Remote Control
- 6 channel MPO for the most advanced fitting options and more individual adjustments

Signia app ready

T-Coil



Signia app and miniPocket
To enhance its discretion, Intuis 3 RIC312 is remotely controllable via the Signia app and the miniPocket.

Intuis 3 RIC312

The remotely controllable Intuis 3 RIC312 is an attractive and convenient solution for mild to severe hearing loss. Thanks to its automatic directionality and a variety of noise reduction features, it offers a reliable hearing experience and improved speech intelligibility even in noisy environments. TeleCare lets your customers enjoy live remote support and fast, smooth trials.

	RIC			
Model				
Receiver type	miniReceiver 2.0 S	miniReceiver 2.0 M	miniReceiver 2.0 P	miniReceiver 2.0 HP
				
Performance				
Matrix (2 ccm)	45/108 dB SPL	60/119 dB SPL	70/124 dB SPL	75/130 dB SPL
Battery				
Size	312	312	312	312
Power consumption	0.9 mA	1.0 mA	1.0 mA	1.1 mA

Proven technology, outstanding comfort

Key features

- For all levels of hearing loss
- Ergonomic design for a high level of wearing comfort
- Remotely controllable for a maximum of convenience and discretion
- TeleCare
- Remote Control
- 6 channel MPO for the most advanced fitting options and more individual adjustments

Signia app ready

T-Coil (except Intuis 3 S)



Intuis 3 BTEs

All Signia's Intuis 3 BTE models offer proven technology with 6 channel MPO for the most advanced fitting options and for more individual adjustments. Suitable for all levels of hearing loss, the easy to handle Intuis 3 BTEs are equipped with automatic directional microphones for a precise hearing experience, offer remote control options, and are TeleCare ready for smoother, faster trials and higher customer satisfaction.

BTE

Model	Intuis 3 S	Intuis 3 M	Intuis 3 P	Intuis 3 SP
				
Performance				
Matrix (2 ccm)	LifeTube 45/124 dB SPL Earhook 55/124 dB SPL	LifeTube 53/125 dB SPL Earhook 60/130 dB SPL	Earhook 70/134 dB SPL	Earhook 80/140 dB SPL
Battery				
Size	312	13	13	675
Power consumption	0.9 mA	1.0 mA	1.4 mA	2.4 mA

The perfect match

Key features

- For mild to severe hearing loss
- 6 channel MPO for the most advanced fitting options and more individual adjustments
- Remotely controllable for a maximum of discretion and convenience
- TeleCare
- Remote Control

Signia app ready

Optional T-Coil



Intuis 3 instant fit ITEs and custom ITEs

With Intuis 3 Instant Fit ITEs and tailor-made ITEs, you can offer your customers the full range of in-the-ear models: No matter if they wish for a solution that fits right away or prefer a customized model, Intuis 3 is the perfect match. Equipped with proven technology, all models offer a precise hearing experience, are remotely controllable for a maximum of discretion and TeleCare ready for successful trials.

Instant Fit ITE / Custom ITE

Model	Intuis 3 Click CIC, Click ITC	Intuis 3 CIC, ITC, ITE
		
Performance		
Matrix (2 ccm)	Click CIC 50/113 dB SPL Click ITC 55/116 dB SPL	CIC 65/124 dB SPL [1] ITC 65/124 dB SPL [1] ITE 65/124 dB SPL [1]
Battery		
Size	Click CIC 10, Click ITC 312	CIC 10, ITC 312, ITE 13
Power consumption	1.1 mA	1.1 mA [1]

(1) Max. values

Basic

Help the world hear better

Basic

The Basic segment of hearing aids offers fitting opportunities that enable you to bring the promise of better hearing to people who would otherwise have no hope of a solution. Easy to learn and easy to fit, even in mobile situations, you can take better hearing to where it's needed – even to outlying areas. This range of BTE and Click ITE devices is suitable for all levels of hearing loss, sizes of ear canal and individual preference.

BTE


Prompt S


Prompt P


Prompt SP


Run S


Run P


RunSP


Fun P


Fun SP


Fast P

Instant Fit ITE


Prompt Click CIC, Click ITC


Run Click CIC


Run Click ITC

ITE


Prompt CIC, ITC, ITE

Make your customers the center of attention

Key features

- A full range of BTE, Click ITE and ITE models
- For all levels of hearing loss
- Remotely controllable for a maximum of convenience and discretion
- TeleCare
- Remote Control

Signia app ready

T-Coil (except Prompt S)

Basic

The Prompt family

Signia's Prompt family offers outstanding wearing comfort and ease of use for every form of hearing loss. All models offer remote control options via the Signia app and the miniPocket™ and are a convenient solution for almost all your customers.



BTE

Model	Prompt S	Prompt P	Prompt SP
			
Performance			
Matrix (2 ccm)	LifeTube 45/124 dB SPL Earhook 55/124 dB SPL	Earhook 70/134 dB SPL	Earhook 80/140 dB SPL
Battery			
Size	312	13	675
Power consumption	0.9 mA	1.4 mA	2.4 mA

Instant Fit ITE / Custom ITE

Model	Prompt Click CIC, Click ITC	Prompt CIC, ITC, ITE
		
Performance		
Matrix (2 ccm)	Click CIC 50/113 dB SPL Click ITC 55/116 dB SPL	CIC 65/124 dB SPL [1] ITC 65/124 dB SPL [1] ITE 65/124 dB SPL [1]
Battery		
Size	Click CIC 10, Click ITC 312	CIC 10, ITC 312/10, ITE 13/312
Power consumption	1.1 mA	1.1 mA [1]

(1) Max. values

Take hearing to a new level

Key features

- Quick and easy fittings everywhere via the MobileFitting app
- For mild to profound hearing loss
- Ergonomic shape for a high level of wearing comfort
- TeleCare
- Remote Control

Signia app ready

Basic

Run

Easy to fit in just a few minutes with the convenient MobileFitting app, Run from Signia is ideal for on-the-spot fittings. Incredibly comfortable thanks to its ergonomic shape, the convenient and easy-to-use Run delivers an outstanding hearing experience.



BTE

Model	Run S	Run P	Run SP
			
Performance			
Matrix (2 ccm)	LifeTube 45/124 dB SPL Earhook 55/124 dB SPL	Earhook 70/132 dB SPL	Earhook 79/140 dB SPL
Battery			
Size	312	13	675
Power consumption	0.9 mA	1.3 mA	2.2 mA

Instant Fit ITE

Model	Run Click CIC	Run Click ITC
		
Performance		
Matrix (2 ccm)	50/113 dB SPL	55/116 dB SPL
Battery		
Size	10	312
Power consumption	1.1 mA	1.1 mA

Provide an outstanding hearing experience

Key features

- A high level of amplification for severe to profound hearing loss
- Ergonomic shape for a high level of wearing comfort
- Simple fittings everywhere via the rocker switch

Basic

Fun

Equipped with Signia's proven hearing aid technology, Fun covers all essential customer requirements. It offers superb sound quality for severe to profound hearing loss, excellent wearing comfort and is simply fitted by pushing the handy rocker switch.



BTE

Model	Fun P	Fun SP
		
Performance		
Matrix (2 ccm)	Earhook 70/132 dB SPL	Earhook 79/140 dBSPL
Battery		
Size	13	675
Power consumption	1.2 mA	2.2 mA

Deliver a powerful listening experience

Key features

- Simple fittings everywhere via the rocker switch
- For mild to severe hearing loss
- Easy handling

Basic

Fast

Fast BTE hearing aids offer a high level of performance and deliver just the right amount of amplification for customers with severe hearing loss. They are easily adjusted via the rocker switch for simplified fitting. Sleek and comfortable with a range of features, Signia's Fast ensures an outstanding hearing experience suitable for every customer.



Model	BTE
	Fast P
Performance	
Matrix (2 ccm)	Earhook 62/130 dB SPL
Battery	
Size	13
Power consumption	1.1 mA



Basic segment feature overview

Available in a wide range of fittings and suitable for all levels of hearing loss from mild to profound, the Basic range offers a good hearing experience for every customer.

Features	Prompt						
	S	P	SP	Instant Fit ITEs	CIC	ITC	ITE
Matrix	LifeTube 45/124 dB Earhook 55/124 dB	70/134 dB	80/140 dB	Click CIC 50/113 dB Click ITC 55/116 dB	65/124 dB [1]	65/124 dB [1]	65/124 dB [1]
Channels/Controls	8/4	8/4	8/4	8/4	8/4	8/4	8/4
Programs	4	4	4	Click CIC: 4 [2] Click ITC: 4	4	4	4
Directional microphones	●	●	●	-	-	-	-
T-Coil	-	●	●	-	-	○	○
Audio input	-	○	○	-	-	-	-
Noise reduction	situation based	situation based	situation based	situation based	situation based	situation based	situation based
SoundBalance™	-	-	-	-	-	-	-
Feedback cancellation	●	●	●	●	●	●	●

Features	Run			Fun		Fast
	S	P	SP	Instant Fit ITEs	P	SP	P
Matrix	LifeTube 45/124 dB SPL Earhook 55/124 dB SPL	70/132 dB SPL	79/140 dB SPL	Click CIC 50/113 dB SPL Click ITC 55/116 dB SPL	70/132 dB SPL	79/140 dB SPL	62/130 dB SPL
Channels	8	8	8	8	6	6	4
Programs	6 preconfigured sound profiles	6 preconfigured sound profiles	6 preconfigured sound profiles	6 preconfigured sound profiles	4 preconfigured sound profiles	4 preconfigured sound profiles	3 preconfigured sound profiles
Directional microphones	●	●	-	-	-	-	-
T-Coil	-	-	-	-	-	-	-
Audio input	-	-	-	-	-	-	-
Noise reduction	situation based	situation based	situation based	situation based	always on	always on	always on
SoundBalance™	-	-	-	-	-	-	-
Feedback cancellation	●	●	●	●	●	●	●

(1) Max. values
(2) 1 program + 4 accessible via app

Color overview SLIM-RIC, RICs and BTEs

Products	Colors												
	Beige	Granite	Grey	Brown	Black	Silver	Pearl White	Sandy Brown	Dark Champagne	Dark Granite	Deep Brown	White	Rose Gold
Styletto IX					●							●	
Motion Charge&Go IX	●				●	●	●	●	●		●		●
Pure Charge&Go IX family	●				●	●	●	●	●		●		●
Styletto AX													
Pure Charge&Go AX / T AX	●				●	●	●	●	●		●		●
Pure 312 AX	●				●	●	●	●	●		●		●
Motion Charge&Go X / P X / SP X	●	●	●		●	●	●	●	●		●		●
Pure 312 X	●	●	●		●	●	●	●	●		●		●
Sirion Connect	●				●								
Orion Charge&Go	●	●	●	●	●								
Intuis 4	●				●	●		●	●				
Intuis 3 RIC 312	●		●	●	●	●							
Intuis 3	●		●	●	●	●				●			
Prompt	●		●								●		
Run	●		●								●		
Fun	●		●								●		
Fast	●												

Products	Colors																
	Dark Brown	Black/ Silver	Cosmic Blue/ Rose Gold	Snow White/ Rose Gold	Black/ Graphite	Black/ Fine Gold	Cosmic Blue/ Silver	Snow White/ Silver	Snow white/ Fine Gold	Graphite	Fine Gold	Cosmic blue/ Fine Gold	Turquoise Blue	Galactic Blue	Pearl Pink	Sporty Red	Color conversion kit
Styletto IX		●	●	●	●			●									
Motion Charge&Go IX										●	●		●	●	●	●	●
Pure Charge&Go IX family										●	●						●
Styletto AX		●	●	●	●	●	●	●	●			●					
Pure Charge&Go AX / T AX										●	●						●
Pure 312 AX										●	●						●
Motion Charge&Go X / P X / SP X																	●
Pure 312 X																	●
Sirion Connect	●																
Orion Charge&Go																	
Intuis 4										●							
Intuis 3 RIC 312																	
Intuis 3																	
Prompt																	
Run																	
Fun																	
Fast																	

Color overview Instant Fit ITEs

Description	Active Pro IX	Silk IX	Silk X	Intuis 3	Prompt	Run
Blue/Mocha		●	●			
Red/Mocha		●	●			
Blue/Black		●	●			
Red/Black		●	●			
Mocha/Beige				●	●	●
Tan/Beige				●	●	●
Black	●					
Snow White	●					
Rose Gold	●					
Beige	●					
Brown	●					
Mocha	●					

Color overview Custom ITEs


Description	Insio Charge-&Go CIC IX	Insio ITC/ITE IX	Insio IIC/CIC IX	Insio C&G AX ITC/ITE	Intuis 3	Prompt
Beige		●		●	●	●
Tan		●		●	●	●
Mocha	●	●	●	●	●	●
Brown		●		●	●	●
Dark Brown	●	●	●	●	●	●
Black	●	●	●	●		

Additional colors may be available upon request.


EarWear 3.0

Discreet. Secure. Comfortable.

miniReceiver 3.0

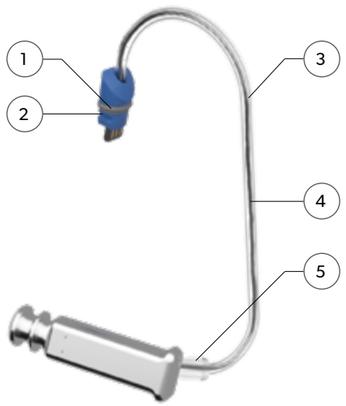

1) Push-/pull system – no tool required
2) Smart Receiver Detection prevents misuse
3) Unique curvature ensues high retention and sits close to head for cosmetic appeal
4) Five lengths 0-4
5) Reinforced cable outlet for high durability
6) Enhanced wax protection
7) Eartip 3.0 and Sleeve 3.0 compatible with miniReceiver 3.0 and ThinTube 3.0
8) WaxGuard 3.0 NanoCare
9) Interface allows easy Eartip 3.0 and Sleeve 3.0 exchange while offering high safety
10) Small interface enables comfortable fit in small ear canals
11) Hypoallergenic Titanium Shield

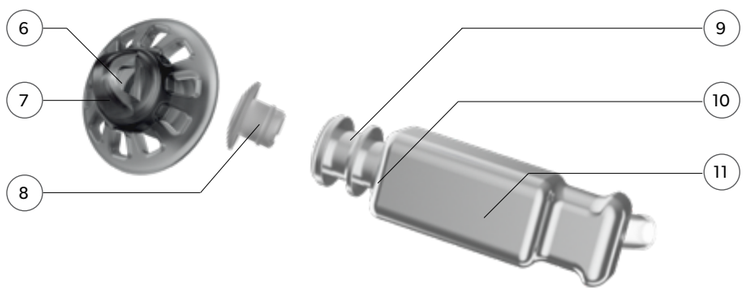

miniReceiver 3.0 available in four power levels: S, M, P, and HP

	S		M		P		HP	
	2ccm coupler	Ear simulator	2ccm coupler	Ear simulator	2ccm coupler	Ear simulator	2ccm coupler	Ear simulator
OSPL 90 (Peak)	110 dB SPL	120 dB SPL	119 dB SPL	129 dB SPL	122 dB SPL	131 dB SPL	131 dB SPL	138 dB SPL
138 dB SPL	46 dB	57 dB	60 dB	70 dB	65 dB	75 dB	75 dB	83 dB


EarWear 3.0

ThinTube 3.0

A choice of two diameters provides a cosmetically appealing solution for mild to severe hearing loss

ThinTube 3.0
0.9 mm

ThinTube 3.0 P
1.4 mm

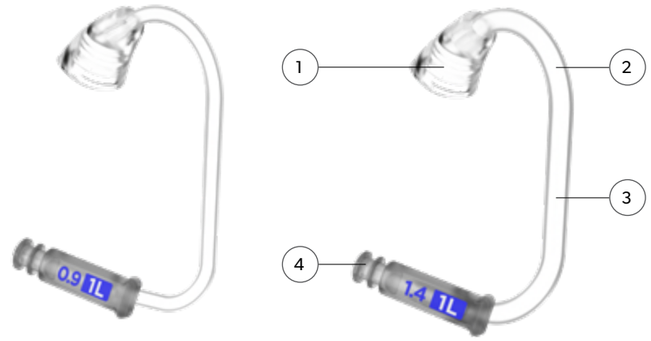

1) One-piece ThinTube 3.0
2) Five lengths 0-4
3) 0.9 mm for mild to moderate and 1.4 mm for severe losses
4) Interface allows easy Eartip 3.0 and Sleeve 3.0 exchange while offering high safety

Eartips 3.0 and Sleeves 3.0


This versatile range is suitable for the vast majority of ear anatomies and hearing needs

Soft eartip and sleeve noses allow for deep insertion with high comfort – also into small ear canals

The secure and proven design makes them much easier for consumers to exchange.

The Power Sleeve 3.0 expands the acoustically highly reliable sleeve to also offer a fully closed option.


All EarWear 3.0 items can be used together.

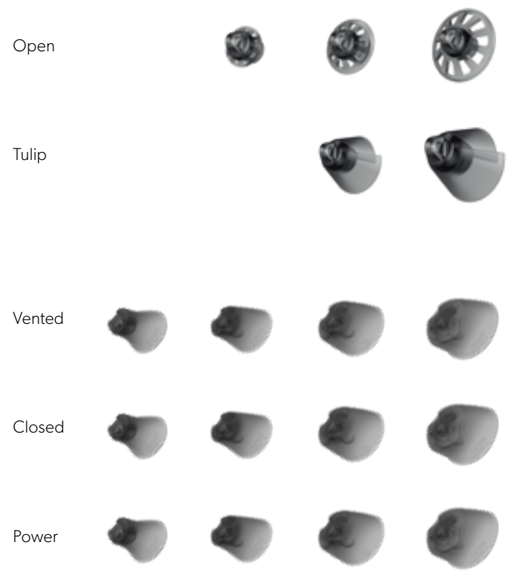

EarWear 3.0

Earmolds 3.0

Compatible with miniReceiver 3.0 and ThinTube 3.0

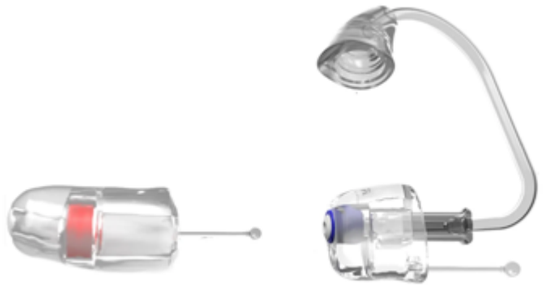

Earmold 3.0 is available in a short and a long version.

The long Earmold 3.0 comes with:
- WaxGuard 3.0 NanoCare
- QuickGuard Wax Guard

Concha Lock Sleeve 3.0

Optional Concha Lock Sleeve 3.0 is available in three sizes



Concha Lock Sleeve 3.0 S – fits miniReceiver 3.0 S and both ThinTube 3.0 and ThinTube 3.0 P

Concha Lock Sleeve 3.0 M – fits miniReceiver 3.0 M

Concha Lock Sleeve 3.0 P – fits miniReceiver 3.0 P

EarWear 3.0

Tools and starter kits

Starter kits are available for: Eartips 3.0, ThinTube 3.0 and miniReceiver 3.0.

Each kit includes specific items relevant to that type of EarWear 3.0, such as:
- Measurement gauge to ensure a perfect fit for miniReceivers 3.0 and ThinTubes 3.0
- Matchbook containing 8 Wax Guard 3.0 NanoCare
- 10-pack of new one-size-fits-all Concha Locks 3.0
- Optional tool for exchanging eartips
- Tool to press out the receiver from a custom eartip

Instruction manuals on how to attach and exchange each component are also included.

EarWear compatibility

	ThinTube 3.0	miniReceiver 3.0	Eartips 3.0 & Sleeves 3.0	Earmold 3.0	miniReceiver 2.0	Click Domes & Click Sleeves	Click Mold
Pure IX series	–	●	●	●	–	–	–
Styletto IX series	–	●	●	●	–	–	–
Active Pro IX series	–	–	●	–	–	–	–
Motion Charge&Go IX	●***	–	●	●	–	–	–
Silk IX series	–	–	●*	●**	–	–	–
Styletto AX series	–	●	●	●	–	–	–
Pure AX series	–	●	●	●	–	–	–
Motion X series	●***	–	●	●	–	–	–
Pure 312 X	–	–	–	–	●	●	●
Silk X series	–	–	–	–	–	●	–
Sirion Connect Series	●†	–	●	●	–	–	–
Orion series	–	–	–	–	●	●	●
Intuis 4 series	●	–	●	●	–	–	–

* Silk Charge&Go IX compatible with Sleeve 3.0 only
** Custom earmolds for Silk IX series not interchangeable with standard Earmold 3.0
*** Motion SP X & IX is not compatible with ThinTube
† Sirion BTE SP is not compatible with ThinTube 3.0

A world of charging possibilities

Portable chargers



Styletto Charger

Pocket-size charging case with integrated power bank that can fully charge a pair of hearing aids up to five times (Styletto IX Charger) or up to three times (Styletto AX charger) before requiring a recharge. Can be charged wirelessly.

Product-specific versions available for Styletto IX, Styletto AX



Pure Portable Charger

Pocket-size charging case with integrated power bank that can charge a pair of hearing aids up to three times before requiring a recharge.

Compatible with all three members of the Pure Charge&Go IX family as well as Pure Charge&Go AX / T AX



Silk Charger

Pocket-size charging case with integrated power bank that can provide up to 3 full charges before requiring a recharge. Can be charged wirelessly.

Use with Silk Charge&Go IX



Active IX Charger

Pocket-size charging case with integrated power bank that can provide up to 3 full charges before requiring a recharge. Can be charged wirelessly.



Portable Charger Insio CIC

A pocket-sized charging case with integrated power bank that can provide up to three full charges for a pair of hearing aids. Can be charged wirelessly.

Use with: Insio Charge&Go CIC IX

Desktop chargers



Dry&Clean Chargers

Not just a charger, it also cleans and dries the hearing aids to ensure peak performance. Achieve a complete UVC cleaning cycle in just 15 minutes.

Product-specific versions available for Styletto IX, Pure Charge&Go IX family, Styletto AX, Pure Charge&Go AX / T AX, Motion Charge&Go P X, Motion Charge&Go SP X, Motion Charge&Go P IX and Motion Charge&Go SP IX



Desktop Chargers

Small desktop chargers with LEDs to indicate charging status.

Product-specific versions available for Pure Charge&Go IX family, Pure Charge&Go AX / T AX, Motion Charge&Go P X, Motion Charge&Go SP X, Motion Charge&Go P IX and Motion Charge&Go SP IX



Desktop Charger Insio AX

Convenient contactless charging means no specific charging position is required - wearers intuitively place their Insio C&G AX hearing aids in the charger.

Use with Insio Charge&Go AX



Multi Charger

One charger conveniently compatible with the following hearing aids:

Pure Charge&Go IX family, Pure Charge&Go AX/ T AX, Motion Charge&Go P X, Motion Charge&Go SP X, Orion Charge&Go, Motion Charge&Go M IX, Motion Charge&Go P IX, and Motion Charge&Go SP IX



Inductive Charger II

Also includes a drying function.

Use with Motion Charge&Go X



Inductive Charger E

For Motion Charge&Go X in performance levels 2X and 1X

Accessory overview

Streamers



The TV Sound streamer provides seamless connectivity, delivering high quality stereo sound for immersive audio. Effortless pairing, intuitive LED indicators, auto-connect functionality, and compatibility with all connectivity-enabled Signia hearing aids.



The StreamLine Mic transforms Signia devices with Bluetooth into more than state-of-the-art hearing aids. They also become hands-free mobile headsets and wireless stereo headphones that stream highest quality audio into both ears from any Bluetooth device. Wearers can even adjust the streaming volume individually via the Signia app.

Remote Control



Small enough to carry on a keyring, the miniPocket is the ideal solution for anyone who wants discreet, convenient hearing aid control without the need for a large remote control or smartphone.

CROS/BiCROS



For people with unaidable hearing loss on one side. Signia CROS devices receive sound on that side and transmit it to a hearing aid on the other side.

App, accessory, and charger compatibility

This overview shows you which apps and accessories work with our hearing aids. Depending on the solution your customers choose, they can even enhance their hearing experience and convenience with this well-matching portfolio.

Hearing aid	Signia app	CROS/ BiCROS	Streamers and transmitters		Remote Controls	Portable Charger	Inductive Charger II	Desktop Charger	Dry&Clean Charger	Multi Charger
		CROS	TV Sound	StreamLine Mic	miniPocket					
Styletto IX	●	●[A]	●	●	●	●	-	-	●	-
Pure Charge&Go IX family	●	●[B]	●	●[F]	●	●	-	●	●	●
Active Pro IX	●	●[C]	●	●	●	●	-	-	-	-
Motion Charge&Go IX	●	●[D]	●	●	●	-	-	●	●	●
Silk Charge&Go IX	●	●[E]	-	-	●	●	-	-	-	-
Insio IX family	●	●[E]	-	-	●	●[G]	-	-	-	-
Styletto AX	●	●[H]	●	●	●	●	-	-	●	-
Insio Charge&Go AX	●	●[I,J]	●	●	●	-	-	●	-	-
Pure Charge&Go AX / T AX	●	●[I,J]	●	●	●	●	-	●	●	●
Pure 312 AX	●	●[I,J]	●	●	●	-	-	-	-	-
Motion Charge&Go X	●	●[K]	●	●	●	-	●	-	-	-
Motion Charge&Go P X	●	●[K]	●	●	●	-	-	●	●	●
Motion Charge&Go SP X	●	●[K]	●	●	●	-	-	●	●	●
Silk X	●	●[L]	-	-	●	-	-	-	-	-
Pure 312 X	●	●[K]	●	●	●	-	-	-	-	-
Sirion Connect	●	-	●	●	●	-	-	-	-	●
Orion Charge&Go	●	-	-	-	●	-	-	-	-	●
Intuis 4	●	●[I,J,M]	●[M]	●[M]	●	-	-	-	-	-
Intuis 3 RIC312	●	-	-	-	●	-	-	-	-	-
Intuis 3	●	-	-	-	●	-	-	-	-	-
Prompt	●	-	-	-	●	-	-	-	-	-
Run	●	-	-	-	-	-	-	-	-	-

A) Compatible with CROS Styletto IX
B) Pure Charge&Go IX / T IX are compatible with CROS Pure Charge&Go IX; Pure Charge&Go BCT IX is compatible with CROS Pure Charge&Go BCT IX
C) Compatible with CROS Active IX
D) Compatible with CROS Pure Charge&Go IX
E) Compatible with CROS Silk Charge&Go IX
F) Not compatible with Pure Charge&Go BCT IX
G) Only for Insio Charge&Go CIC IX
H) Compatible with CROS Styletto AX
I) Compatible with CROS Pure Charge&Go AX
J) Compatible with CROS Pure 312 AX
K) Compatible with CROS Pure 312 X
L) Compatible with CROS Silk X
M) Excluding PL 4.0
N) Compatible with Bluetooth-enabled Insio Nx ITC/ITE

All your clients' needs, all together



Signia app
Greater convenience

Simple and convenient, the Signia app empowers wearers with access to remote control adjustments, audio streaming, and the Signia Assistant. And because it includes TeleCare, it integrates seamlessly into your fitting process, making it easier for your customers to stay in touch with you when they are away from your premises. Suitable for use with all Signia hearing aids, its smart user interface only shows the user functions applicable to their own devices.



Download on the App Store | GET IT ON Google Play

Compatible with Android and iOS.
The app is free of charge.

TeleCare

Stay in closer contact with your customers:
- Real-time contact via text message, voice call or video call
- Intuitive monitoring of your customers' daily satisfaction
- Set virtual appointments for the precise remote adjustment of the hearing aids

Remote Control

Personalize your customers' hearing experience further:
- Hearing program and volume control
- Sound Balance Adjustment
- 360° Spatial Configurator for individual adjustment of directionality
- Continuous monitoring of the hearing aids' connection and battery status

Connectivity & Streaming

- Easy management of audio streaming accessories (e.g. StreamLine TV, StreamLine Mic) as well as CROS devices
- Volume level control while streaming
- Adaptive streaming volume

TeleCare | Remote Control | Connectivity & Streaming | Signia Assistant

Signia Assistant
The smart new way to better address your clients' needs


Imagine being able to help your clients to Be Brilliant every single step of their hearing journey outside your office? Now you can – with the revolutionary Signia Assistant.


The world's first live deep neural network for hearing aids

Easy to use, the Signia Assistant offers small sound adjustments when needed and how-to videos to answer handling questions. With its live artificial intelligence deep neural network, the Signia Assistant gives your clients the best possible solution when issues arise in between appointments. All they need is their Signia IX or AX hearing aids and the Signia app.

See every detail with Connexx

With the Connexx fitting software, you'll find all the changes made by the Signia Assistant, enabling you to understand each patient's key challenges. It's the smart solution for more efficient consultations and fittings - and more satisfied customers.

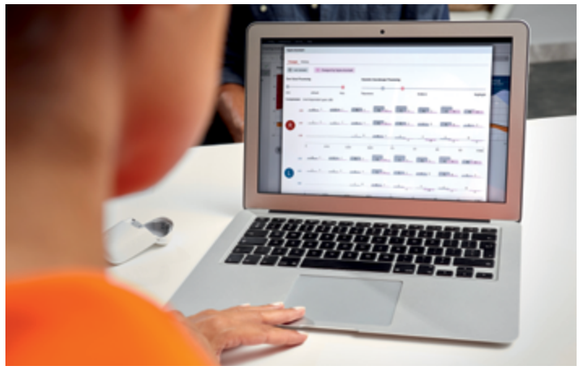

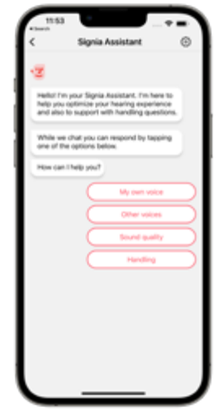

With the Signia Assistant you can:

- Give your clients the assurance of 24/7 support
- Increase your client satisfaction and conversion rate with the power to address hearing aid wearers' needs more accurately
- Increase the quality of follow up appointments through a new level of real-world usage insights
- Fine-tune with ease by seeing exactly how your clients choose to personalize their hearing experience

Successful fittings, satisfied clients



Connexx Fitting Software



Unity™ 4
mobile fitting unit system

Signia offers you a broad range of counseling and fitting support – from programming software and tools to mobile fitting units to hearing aids for demonstration and advisory tools – up to comprehensive remote hearing care. This enhances your workflow and makes trials and fittings smoother and more pleasant for your customers.

Fitting Integrated Xperience in Connexx

Our Connexx fitting software puts you in the driver's seat with real-world insights to optimally support each individual patient.

First Fit in Connexx configures Signia IX hearing aids to allow the wearer to hear what is important to him or her. The First Fit settings have been carefully optimized to provide a good starting position for each wearer.

The Dynamic Soundscape Processing (DSP) control in Connexx provides a slider to adjust to your clients' individual preferences. This slider is now more powerful than ever thanks to our new Augmented Focus™ technology.

Unity 4 hearing assessment and fitting unit

The Unity 4 fitting unit boasts a small physical footprint, showcasing our revolutionary stackable design. Embrace the perfect blend of functionality and space-saving elegance.

Elevate the client experience, enhance workflow, and confidently deliver superior care with our hearing assessment and fitting unit.

Programming Interface Compatibility

Hearing Aid	Noahlink Wireless	Connexx Air	I3 Adapter	Flex Connector	Size 312 Adapter	Size 13 Adapter	Programming Shoe Size 13	Programming Shoe Size 675	MobileFitting app	Fit via rocker switch
Styletto IX	●									
Active Pro IX / Active IX	●									
Motion Charge&Go IX	●									
Pure Charge&Go IX family	●									
Silk Charge&Go IX		●								
Insio IX family		● [1]	● [2]	● [3]						
Styletto AX	●									
Insio Charge&Go AX	●									
Pure Charge&Go AX / T AX	●									
Pure 312 AX	●									
Motion Charge&Go X / P / SP X	●									
Silk X		●		●						
Pure 312 X	●									
Sirion Connect	●									
Orion Charge&Go	●									
Intuis 4 BTEs	●									
Intuis 3 RIC312					●					
Intuis 3 S					●					
Intuis 3 M						●				
Intuis 3 P							●			
Intuis 3 SP								●		
Intuis 3 Click CIC / ITC				●						
Intuis 3 CIC, ITC, ITE				●						
Prompt S					●					
Prompt P							●			
Prompt SP								●		
Prompt Click CIC / ITC				●						
Prompt CIC / ITC / ITE				●						
Run S / P / SP (BTEs)									●	
Run Click CIC / ITC									●	
Fun (all)										●
Fast										●

1 Versions with e2e wireless technology only
2 Insio IIC/CIC IX with I3 faceplate, required for non-wireless versions and compatible with e2e wireless versions
3 Insio CIC IX with T5 faceplate and Insio ITC/ITE IX, required for non-wireless versions and compatible with e2e wireless versions